



VidimScan® Universal Gel ID Card-Cassette Reader IVD approved



Features :

VidimScan® is universal: processes all Gel ID Card –Cassettes on the market.

VidimScan® will give a interpretation of the blood group results.

VidimScan® will read, scan, archive , digitalize, print ,give full color image.

VidimScan® has unique security features = Patient barcode info is linked to the barcode of the Gel Id Card-Cassette.





Ortho, Biorad, Grifols, Diamed, Tulip, Invitrogen, ...

VidimScan® Universal Gel Card-Cassette Reader can be used to process all gel cassettes for blood group serology tests on the market today . Our frames are compatible with all gel cassettes on the market . Depending of the chosen type of cassette you can put 9 or up to 12 cassettes in one frame

The IVD approved software of VidimScan® contains the parameters of the following tests:

ABO/Rh , ABD-confirmation , ABO/Rh for newborns , ABO/Rh for patients , ABD-confirmation for patients , ABO/Rh for newborns DVI+ , ABO/Rh for newborns DVI- , ABO Rh for Donors, ABD- confirmation for donors , ABO/RhD , ABO/D , ABO/D+DAT, ABO/D+ reverse grouping, Rh-subgroupes+K, Rh-Subgroups + Cw +K, RhD + Phenotype, Anti-A1 Absorbed, Anti-H, ID-Partial RhD typing, Coombs Anti-IgG, ID-Antigen Profile I, ID-Antigen Profile II, ID-Antigen Profile III, Anti-D, Anti-DVI neg, Anti-Cw, Anti-K, Anti-K/k, Anti-M, Anti-N, Anti-M/N , Anti-P1, Anti-Le a, Anti –Le b, Anti-P1, Anti-Le a / Le b, Anti-Jk a, Anti-Jk b, Anti- Jk a / Jk b , Anti-Kp a, Anti-Kp b, Anti-Kp a/ Kp b , Anti-Lu a, Anti-Lu b, Anti- Lu a/Lu b, Test of Antigen S, Test of Antigen s , Test of Antigen S and s , Test of Antigen Fy a, Test of Antigen Fy b, Test of Antigen Fy a and Fy b, ID-Anti-I, ID-Anti-i, Reverse grouping with Antibody Screening, NaCl.Enzyme Test and Cold Agglutinins, LISS / Coombs, Coombs Anti-IgG, DC Screening I, DC Screening II, DAT IgG-Dilution, DT IgG1/IgG3, Screen, LISS :Coombs +Enzyme test, LISS/Coombs, Coombs Anti-IgG, NaCl .Enzyme Test and Cold Agglutinins, Type + screen, Complete Crossmatch, LISS/Coombs, Coombs Anti-IgG, NaCL ;Enzyme Test and Cold Agglutinins, ABD-Confirmation, AB- confirmation for Patients and for Donors....

Update :

Vidimsoft software is constantly updated with new tests and additional tests and information from Gel ID Card –Cassette producers.

Patent Pending

Vidimsoft software is protected by patents



VidimScan® can be used as :

a standalone equipment for manually operating labs

or to process blood from patient after blood transfusion (mixed results)

or to process baby blood tubes which do not fit in some automates

or as internal quality control, with traceability

or as emergency backup for fully automated labs

or as a quick alternative when the automate is in full process




Unique security features

Unmatched today on the market are Vidimscan security features .
Security and avoiding mistakes is a top priority in immunohematology .

Unique security combination

Link between barcode Patient & Gel ID Card / Cassette

barcode gel Id card

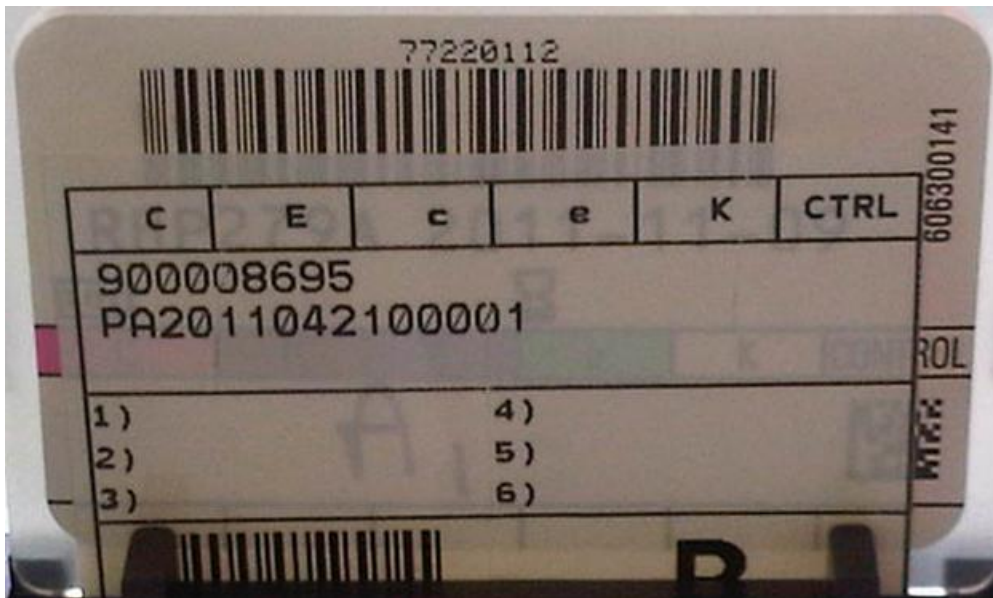


barcode patient information

Patient coordinates/number/blood tube sample / barcode are linked with barcode Gel ID Card-Cassette.

Database for traceability

VidimScan® will give a warning for example when actual input doesn't fit with blood tube ID history database


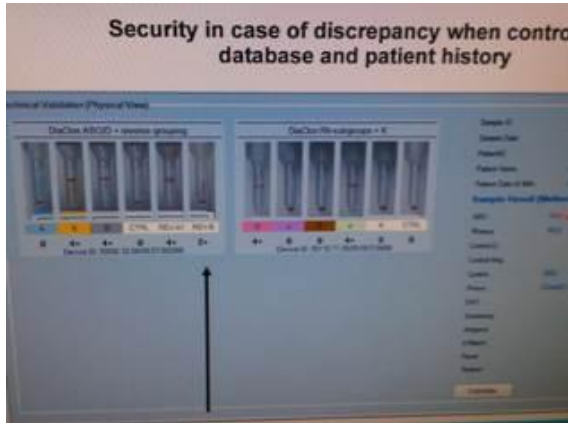


Final important note : Blood group plates

VidimScan® can be combined with VidimPlate® .

This will allow two operators to work on the same patient blood tube , independently from each other with two different techniques.

See VidimPlate® for more information

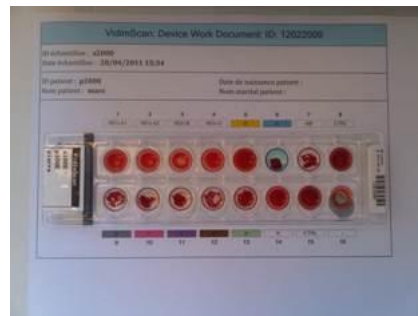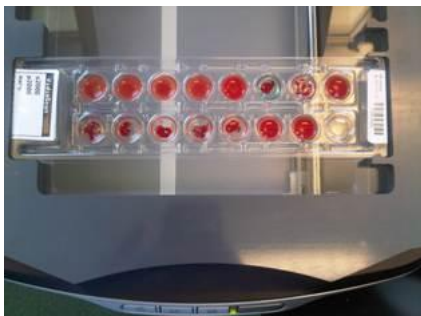

VidimScan® Universal Gel ID Card Reader

Soft-and hardware for digitalization & interpretation & traceability for manual blood group serology test .

Software IVD approved , meets Cofrac 15189, meets ISO15189 and Directive N° 2005/36/CE 7.9.2005 EU.

A stand alone device digitalizing all gel cassettes Ortho, Biorad, Grifols, Diamed, Tulip, ...

Full color image archiving and consulting, and enlargement of the image .

Database for traceability of all performed handlings during test (name, test, date ...)

Reagents tracking

Unique security features , link between patient and test

User friendly

Full patient blood tube history

Data exchange interface with the LIS, LIMS

Personalized record

Built in quality control system

Software can be translated into your own language (for partners)

Tele-maintenance support , through internet access.

24/7 service backup worldwide

Can also be used for processing of plates , or a combination of plates and Gel Id Card-Cassettes

VidimScan® brings the benefits of digitalization and traceability and security to manual blood grouping.

VIDIMSOFT bvba *your security on blood grouping*

Veldstraat 33 B-8850 ARDOOIE Info : +32.51.54.12.00 Mobile : +32.473. 93. 38 .36 Fax : +32.51.54.65.88

Email : marc.desomer@vidimsoft.com Website : www.vidimsoft.com